



17 x 51 in. Long-Format DR Wall Stand

One exposure, one image, no stitching

Model: 17 x 51 Long Format Wall Stand (for DFMT Elite and OTS systems)



This wall stand captures a full 17 x 51 in. image in one exposure, ideal for full spine, scoliosis and long leg studies. For routine exams, use just the top, middle or bottom section as a 17 x 17 in. image. It can be part of a new ceiling-mounted or floor-mounted room, or be added to an X-ray room you already have.

- **No stitching**

One exposure covers the full spine or both legs, with no image matching and no patient movement between exposures.

- **Three fields of view**


17 x 51 in., 17 x 34 in. or 17 x 17 in.

- **Easy to position**

Dual lock release handles and dual remote collimator light controls.

- **Retrofit option**

Add to a new room or an existing X-ray suite.

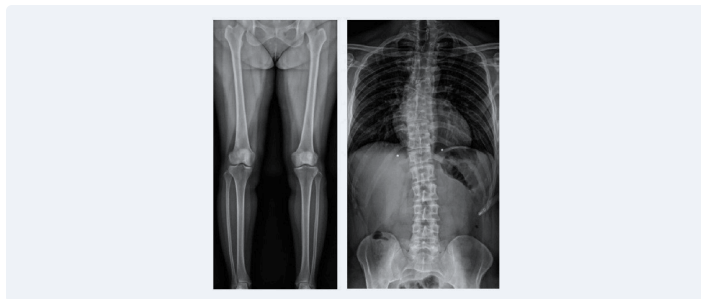

Long leg and full spine images

Specifications

Long-format detector		Environment	
Image sensor	a-Si TFT	Operating	+10 to +35 °C; 30% to 85% humidity, non-condensing
Scintillator	Gd2O2S:Tb (gadolinium oxysulfide)	Storage and transport	-15 to +55 °C; 10% to 90% humidity, non-condensing
Pixel pitch	140 µm		
17 x 51 in. field	Active array 3072 x 9216 pixels		
17 x 34 in. field	Active array 3072 x 6144 pixels		
17 x 17 in. field	Active array 3072 x 3072 pixels		
Grayscale	16-bit		
Spatial resolution	3.5 lp/mm		
Image time, 17 x 51 in.	6 s preview, 12 s final; 5.8 s cycle		
Image time, 17 x 17 in.	3.5 s preview, 5.9 s final; 3.7 s cycle		
X-ray sync	AED (automatic exposure detection)		
Image transfer	Gigabit Ethernet (1000BASE-T)		
Size and weight	1357 x 532 x 30 mm; 18.5 kg		

Sold, installed and serviced by Sharp Imaging Solutions

- **Site planning and installation.** Room review, deinstallation of old equipment, installation and calibration.
- **Canon CXDI DR.** Authorized Canon dealer and service provider for detectors, software and DR upgrades.
- **Connectivity.** DICOM worklist and PACS setup, with training for your technologists.
- **Service and parts.** Preventive maintenance, repairs and parts for multi-vendor X-ray equipment.

Request a quote

561-441-2157

sales@sharpimagingolutions.com
sharpimagingolutions.com

Manufactured by AMRAD Medical, a brand of Summit Industries, LLC. Sharp Imaging Solutions LLC is an independent dealer that sells, installs and services this equipment, and is not the manufacturer.

Specifications are based on the manufacturer's published information and may change without notice. Final configuration, specifications and warranty terms are confirmed in your written quote.

AMRAD Medical, DFMT Elite, OTS and related product names are trademarks of Summit Industries, LLC. Canon, CXDI and related names are trademarks of Canon Inc. Windows is a trademark of Microsoft Corporation. These names are used only to identify the equipment and do not imply endorsement.

Product photos are the property of the equipment manufacturer and are used with permission.

© 2026 Sharp Imaging Solutions LLC. Brochure text and layout by Sharp Imaging Solutions LLC.